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Thank you for choosing this product.
Please read this manual carefully before using it. Please keep it in good condition for future inquiries.

Product Efficiency

- Relieving muscle fatigue and soreness
- Release tension
- Promoting physical and mental relaxation
- Promote blood circulation
- Reduce fatigue and burnout
- Improve the sub-health status of the cervical spine

Parameter Description

- Product name: Massage pillow
- Rated power: 3.7W
- Automatic timing: around 10 mins
- Charging voltage: DC 5.0V
- Net weight: 495g

Usage Method

1. To turn the device on press and hold the button III for 2 seconds.
2. Press the M button to switch between the modes: massage, reverse massage, vibration, heating.
3. Press the button III to change the intensity of the massage.
4. To turn the device off press and hold the button III for 2 seconds.
5. The massager is powered by a lithium battery. If you notice that the unit's output has significantly decreased, charge the massager.
6. When charging, please use the charging cable of this product to check whether the power connection is normal, plug in the DC interface to charge.
7. When charging, the voltage may be unstable, and the product is forbidden during charging.
8. It takes about 2 hours to charge, and the power is unplugged after being fully charged.
9. This product includes accessories only for use with this machine. Other market accessories are not a substitute for this product.


Common problems and handling

| Serial number | Phenomenon | Possible Reason | Approach |
|---------------|----------------------------|---|---|
| 1 | Massager can't start | Exhausted during standby | It can be started normally after charging |
| 2 | The massager stops working | The timing massage program has been completed or the battery power has been used up | Re open the battery switch |


If there is still an abnormality in the above processing, please contact your local dealer for repair service.

Product Configuration


Your U-shaped pillow includes: main unit, charging cable, instruction manual, certificate

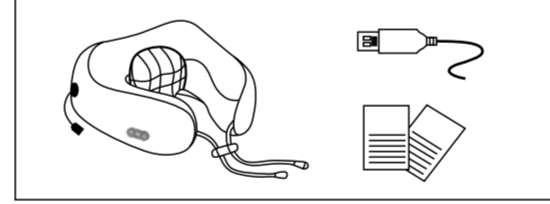Kneading type

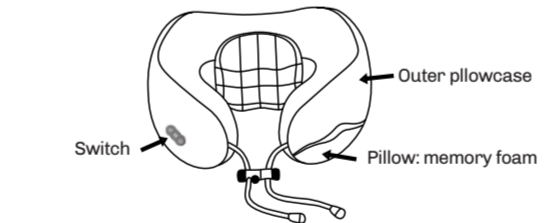
Rope buckle design

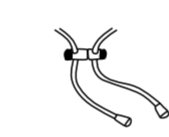
Press the unlock button

Press up and down to loosen

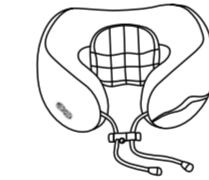One button start
Switch modes
Change intensity

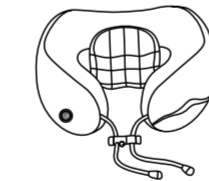One button start

Use The Power Plug And Powercord

1. Please check that the voltage does meet the product specifications.
2. Do not use other similar accessories instead of this product accessories.
3. Do not touch the power supply or plug with wet hands to avoid malfunction or electric shock.
4. Do not charge the product when the voltage is unstable.

Do not use this product when the following conditions occur:

1. The massager is abnormal, and the massage device is damaged or damaged.
2. If there is any serious or severe pain during the massage, stop the massage immediately and seek medical advice.
3. In case of the above situation, if it is confirmed that the product is the cause, please send the product to the repair service center for inspection and maintenance.

Important Safety Instructions

Before using this product, please read all the materials and follow the instructions below:

1. Do not place the product in the bathroom or in a humid place.
2. Do not leave the massager in direct sunlight or high temperature for a long time.
3. If the product is used by children, physically and mentally disabled or disabled, it needs to be accompanied and guided by the guardian.
4. This product has a built-in lithium battery, please dispose of this product in accordance with the relevant environmental protection in the city.

Not Suitable For Use Of This Product

If you have any of the following conditions, please consult your doctor before use.

1. Being treated by a doctor or physically discomforted
2. Patients with acute, sudden patient
3. Drinking users
4. Minors and teenagers, children

Battery Charging Use

1. Under normal use, the life of the rechargeable battery is 2500 times charge and discharge.
2. When the rechargeable battery is not used for a long time, the battery life will be shortened. It is recommended to charge at least 1-2 times a month.